

Thank you for purchasing this product from Gopher.



**UNCONDITIONAL
100% SATISFACTION
GUARANTEED**

If you are not satisfied with any Gopher purchase for any reason at any time, contact us and we will replace the product, credit your account, or refund the purchase price.

Read instructions in their entirety. In the event of missing parts, call **1-800-533-0446**.

Mac[®] Software Manual for FITstep[™] Pro Version 3

Operating System Requirements

Item	Minimum	Recommended
Operating System	Mac [®] OS X [®] 10.9 Mavericks	Mac [®] OS X [®] 10.11 El Capitan
Processor	Intel [®] Core [™] 2 Duo, Core [™] i3, Core [™] i5, Core [™] i7, or Xeon [®] Processor	Same
Memory Capability	2 GB	Same
Free Hard Drive Space	100 MB (including installed program and database)	500 MB (including installed program and database)

Installing the Software

To download the FITstep[™] software, go to www.gophersport.com/software. Click on the Mac Version 3 User Software Download, follow the prompts to install.

Drag the FITstep[™] Pro logo to the Applications folder. As part of the installation process, any data you collected with the first version of the software is converted to work with the new version. It is not necessary to re-enter student names, pedometer number, or classes. Also it is important to note, that the installation process will create a backup file that you can use with your older version of the FITstep[™] software.

Using the Software Menu

The screen pictured on page 2 appears when the software opens. The components in the menu bar are used for file operations, printing, and importing/exporting data.

The screen is divided into two major parts. The left window is hierarchical and breaks down by instructor to grades to classes to students. To make the next level appear, double click on the top level tile or click on the triangle to the left of the name. The right window is the main viewing window; class/student information is entered here, full class lists are shown, and uploaded data can be viewed. To use any of the tabs on the right side of the screen, a class must be selected.

It is possible to utilize the Reports Tab function with more than one instructor, student, class, and/or grade by holding the command key down and clicking the names of the desired items.

Follow these steps to create a new database:

1. File
2. New Database
3. Select the location you wish to store the database (desktop is recommended.)
4. The software will prompt you to name the database.
5. Select Create new FITstep Database.

Once database is located or created, proceed to Reader Setup Instructions.

Reader Setup

The reader allows for data transfer from the FITstep™ Pro pedometer to the FITstep™ software.

1. To set up the reader, plug the reader into the computer using the provided USB cord.
2. The computer will automatically connect to the reader.
3. The FITstep™ software will register the reader as indicated in the lower left hand corner.

TABS

Pedometer Setup

The Pedometer Setup Tab has one purpose: to embed a number in the pedometer. Embed pedometer numbers before using them with students. Here is the suggested procedure for numbering the pedometers if you have a new set:


1. Number the pedometers on the outside and inside of the pedometer covers. This will make it easy for students to identify their pedometer each time they come to class. It is probably easiest to start with pedometer 1 and number from that point forward. However, numbering can be started at any integer.
2. Pedometer may be in a deep sleep mode to conserve battery life. To wake the pedometer, press and hold the reset button until the screen comes on.
3. Click the Pedometer Setup tab. Click on the “write new number” box. Set the new pedometer number to the desired starting number. For example, if you set the new number to 1, that pedometer will be internally assigned to number 1.
4. Select the first pedometer you wish to number, open the pedometer and place it screen side up into the reader.
5. The software will register the pedometer and assign it with the number identified.
6. The software will automatically change the next pedometer number to the next integer. You are now ready to open the next pedometer and place it into the reader.
7. To check whether the correct number has been assigned to the pedometer, press the Mode button on the pedometer until the Step screen is showing. Now hold down the Set button to see the pedometer number. It will appear as P=00001 if you embedded number 1.

NOTE: nothing will show on this screen if the pedometer has not been numbered.

8. To manually change a pedometer number, check the “Write New Number” box and type in a new number. Open the front cover of the pedometer, place it face down into the reader and the new number will be assigned. A pedometer number can be changed using only the Pedometer Setup Tab.

Class Setup

The Class Setup Tab serves as the control center where you will add new instructors, grades, classes, and students. Not all fields will need to be filled, such as the instructor number and description fields. Notice how the data is organized. The top field is the instructor, usually the Physical Education teacher. Next is Grade followed by Class Name, which might be listed as a classroom teacher in an elementary school or as a class number in secondary schools (for example 3rd Period). If desired, type a description of the class with information that might be relevant. Finally, in the last section add the students by first and last name, gender (optional), and pedometer number. If desired, a student ID can also be added.

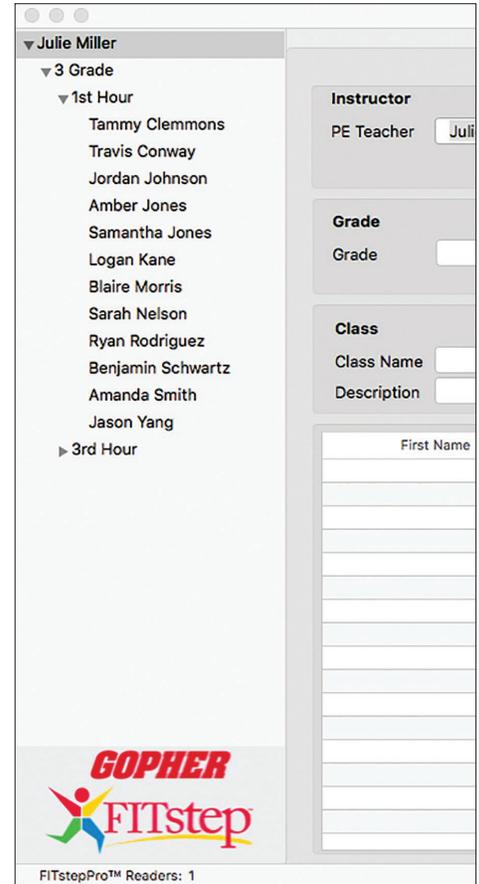

FITstep™ v3.3

Class Setup
Upload Pedometer Data
Pedometer Setup
Reports

Julie Miller

▼ 3 Grade

▼ 1st Hour

- Tammy Clemmons
- Travis Conway
- Jordan Johnson
- Amber Jones
- Samantha Jones
- Logan Kane
- Blaire Morris
- Sarah Nelson
- Ryan Rodriguez
- Benjamin Schwartz
- Amanda Smith
- Jason Yang

► 3rd Hour



FITstepPro™ Readers: 1

Instructor

PE Teacher Instructor No.

Grade

Grade

Class

Class Name

Description

First Name	Last Name	Gender	Pedometer Number	Student Num...
Tammy	Clemmons	Female	11	
Travis	Conway	Male	10	
Jordan	Johnson	Male	2	
Amber	Jones	Female	5	
Samantha	Jones	Female	3	
Logan	Kane	Male	12	
Blaire	Morris	Female	6	
Sarah	Nelson	Female	8	
Ryan	Rodriguez	Male	9	
Benjamin	Schwartz	Male	4	
Amanda	Smith	Female	1	
Jason	Yang	Male	7	

Manual Entry Class Setup

Teachers, students, and classes can be entered into the FITstep™ software manually by typing the information into each specific field.

Import Students into Class

The Import/Export menu offers the opportunity to import students into the FITstep™ software using a spreadsheet.

NOTE: Data can only be uploaded from an Excel 2007 or higher (.xlsx) workbook. Multiple classes can be stored in the same workbook, but each class must be in a separate worksheet. The worksheets do not need to be named, however naming them makes it easier to identify different classes that exist in the same spreadsheet.

The workbook should be set up in columns as follows: First Name, Last Name, Gender, Pedometer Number, and Student Number (see example below).

NOTES:

- Data cannot be imported unless a first and last name is provided.
- Pedometer numbers must be integers greater than zero.
- The software will accept gender identifiers of "Male", "M", "Boy", or "B" for males.
- The software will accept gender identifiers of "Female", "F", "Girl", or "G" for females.
- The student number field is available for adding student IDs or another unique identifier. This field accepts characters and numbers.
- Columns do not need to be in any particular order.

	A	B	C	D	E	F	G	H	I	J	K
1	First Name	Last Name	Gender	Pedometer Number	Student Number						
2	Katie	Albertson	F		1 A487361						
3	Andy	Alito	M		2 C845963						
4	Katelyn	Anderson	F		3 A487361						
5	Carter	Aoki	M		4 C845963						
6	Molly	Brickley	F		5 A487361						
7	Jessie	Collins	M		6 G697435						
8	Matthew	Cranson	M		7 B123456						
9	Monica	Draben	F		8 A487361						
10	Ellie	Elias	F		9 T741596						
11	Alexina	Garza	F		10 G697435						
12	Leo	Gonzalez	M		11 G697435						
13	Nick	Jackson	M		12 B123456						
14	Amy	Marbut	F		13 T741596						
15	Donald	Omerta	M		14 G697435						
16	Aiden	Pryor	M		15 B123456						
17	Emma	Schaeffer	F		16 T741596						
18	John	Smith	M		17 B123456						
19	Isabelle	Wiggins	F		18 T741596						
20	Juan	Wilson	M		19 C845963						
21	Braden	Woods	M		20 C845963						
22											

Export Full Database to Excel

The Import/Export menu also allows you to place all data into an Excel workbook. The created workbook will contain five worksheets for instructors, grades, classes, students, and workouts. This command should not be used for backing up data since once it is exported to Excel, it cannot be imported back into the FITstep™ software. The primary purpose for this command is to do further statistical analysis in Excel.

Upload Pedometer Data

The software works with the FITstep™ Pro data readers. The reader quickly uploads the data (about 1.5 seconds) per student. Multiple readers are useful because a large class can be uploaded rapidly. You must plug readers into a USB port on your computer. If you need to use multiple readers, a USB hub is recommended. The number of readers actively connected to your computer is listed in the bottom left of the screen.

NOTE: Ensure that each student clears his/her pedometer prior to beginning the class activity. This is done by pressing and holding the Yellow Reset button. Pedometers do not automatically clear data from previous uses.

1. Select the class on the left hand side of the screen.
NOTE: Data cannot be transferred from class to class.
2. Select the Upload Pedometer Data tab.
3. Students open the front cover of the pedometer and place it screen side up into the reader.
4. As data is uploaded, the red and green lights on the reader will flash and students' names in the left window will change color and re-sort as shown in the screenshot below. The color of a student's name offers the user the following information:
Black Font: students have not put their pedometer into the reader.
Green Font: data was successfully uploaded.
Red Font: data upload was unsuccessful. This may occur when the student does not leave the pedometer in the reader long enough for the data to upload.
Blue Font: uploading process is in progress.
5. When uploading is complete, only the red light stays solidly lit.

NOTE: The key to successfully uploading data is to leave the pedometer in the reader until the red and green lights stop flashing.



My Physical Activity Report

Our school collects physical activity (PA) data using pedometers. Pedometers accurately monitor the amount of activity students collect in three different settings. First, we collect activity data in physical education classes. When students are comfortable using the pedometers during physical education, they wear them throughout the school day. In the future, those students who demonstrate responsible care of the pedometers may be allowed to take the pedometer home to see how much activity they accumulate daily (24-hour period). Typically, four days of full-day data is gathered for elementary school students and eight days for secondary students. The days do not have to be consecutive so the time period during which the data is collected may be greater than 4 - 8 days.

Student: Amanda Smith
Class: 1st Hour
Grade: 3 Grade
Teacher: Julie Miller
Data Gathered: 2/3/18 to 2/14/18

	Student Total	Student Average	Class Average
Steps:	48	48	17
Activity Time	0:20	0:20	0:06
MVPA Time:	0:13	0:13	0:04
Days Measured:	1		

Steps: The number of steps accumulated during the time period listed above.

Activity Time: The total number of minutes accumulated regardless of how fast or intense the physical activity is.

MVPA: An acronym for Moderate to Vigorous Physical Activity that is defined by various educational agencies as a brisk walk at 3 to 4.5 miles per hour. Students receive more health benefits from MVPA.

Print Reports

The print function is listed under the File menu. It is available when students have been selected and the data is listed on the screen.

Click on the print command to see a preview of your report. If you are satisfied with the report, click on the print icon to print the report. If not, cancel the report and filter the data as preferred.

Another feature of the Report Viewer is that the report may be saved in many different formats, including PDF or Word Document. To do so, click on the triangle to the right of the disc (see arrow) on the command bar of the report and pull down to your desired formatting.

Monthly Summary

Instructor	Grade	Class Name	First Name	Last Name	Ped #	Time Range	Days	Steps	Avg Steps	Activity	Avg Activity	MVPA	Avg MVPA
Julie Miller	3 Grade	1st Hour	Samantha	Jones	3	2018-02-01	1	39	39	0:16	0:16	0:08	0:08
Julie Miller	3 Grade	1st Hour	Benjamin	Schwartz	4	2018-02-01	1	49	49	0:17	0:17	0:14	0:14
Julie Miller	3 Grade	1st Hour	Amanda	Smith	1	2018-02-01	1	48	48	0:20	0:20	0:13	0:13
Julie Miller	3 Grade	1st Hour	Jordan	Johnson	2	2018-02-01	1	56	56	0:23	0:23	0:17	0:17

Reports Prepared by FITstep™, Gopher Sport © 2018

1

Export Report to Excel

A report may be exported to Excel by selecting this option under the File menu. This allows any report chosen to be exported into an Excel spreadsheet. The spreadsheet will be stored in the same folder where the data file is stored. If the file hasn't been moved since startup, it will be stored in the folder, c:\ProgramData\Gopher\FITstep\.

General Report Help

All the reports generally function in a similar fashion. You pick the students to report on from the class tree on the left side of the main window. You can select the instructor, grade, class, or student level. Remember that the Command Key can be used to select a number of distinct items in the tree or the Shift Key to select ranges of data. Next, select a date range for the report. The report can be customized by selecting options that are available on the right hand side of the main window. Finally, click on the column name to sort the data. The first click will sort the data in ascending order, the second click will sort in descending order and the next click will return the data to its original order.

After selecting the students/classes/grades/instructors and report features, print the report by selecting "Print Report" from the File menu. This generates a separate report viewer window. From this window it is possible to print the report to a printer, view formatted pages on screen, or export the report to one of several file formats. To print, simply click on the printer icon and follow the steps in the printer setup menu. Use the navigation arrows to scroll through the report. The file can also be saved to any of the following file formats: Acrobat (PDF) File, CSV (comma delimited), Excel Worksheet, PowerPoint Presentation, Rich Text Format (RTF), TIFF File, Web Archive, Word Document, or XPS Document.

It is possible to "Print from Excel" from the File menu. This creates an Excel 2007 format worksheet that has the data items that are selected for the report. The Excel file that is created is a data dump of items that you have selected in the order you have chosen.

Preferences Menu

The Options Menu allows the user to set some defaults for the software.

Select Start Date. This option is used for selecting a starting date for creating reports and screen analysis of data. This command dictates what the starting date will be when the software is launched. It will always start at one of the options listed below until the setting is changed. The following date range options may be selected:

Today – defaults to today.

Yesterday – defaults to yesterday.

Current Week – defaults to Sunday of the current week.

Current Month – defaults to the first day of the current month.

Current Year – defaults to January 1 of the current year.

Last Week – defaults to Sunday of the previous week.

Last Month – defaults to the first day of the previous month.

Fixed Start Date – defaults to a fixed date that is selected. This day will not change until you adjust it. You can use this feature to set the default date to the start of the current term, if desired.

Custom Start Date – allows data to be viewed from a certain number of days, weeks, months, or years. Select the time period (days, weeks, months, years) to go back and then select the number of time periods you wish to go back. When using this option, you will not necessarily go back to the start of a week, month, or year. When going back by weeks, it will default to the same day of the week as today. When going back by months, it will default to the same day of the month as today. Finally, when going back years, it will default to the same day of the current month in the previous year.

After choosing an option, click OK to verify the selection or click Cancel to abort any changes. The next time the Reports tab is selected, the start date for the reports will default to the option selected.

Show Advanced Data Options. This option can be used to check the amount of time students spent in each of the steps per minute levels. The FITstep™ Pro pedometers show how much time was spent at each SPM levels (80 to 160 steps per minute). By selecting the advanced data option, this data can be reviewed.

NOTE: This feature only works when viewing the Workouts tab in Reports.

Customize Activity Reports. This option works on the Physical Activity report and can be used to customize reports. It allows changes to be made to the Physical Activity report to suit the school preferences. The changes can be reviewed using the preview selection.

The options that can be changed are:

- Edit Report Title
- Edit Intro Paragraph
- Change School Logo
- Add Logo to Report

Help Menu

The Help Menu provides easy access to the FITstep™ Software instructions for reference during use. Additional information can be found by following the Gopher Sport Web Links.

Future updates to the FITstep™ Software will be available for download at www.gophersport.com/software.